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Important note regarding use

Please read these instructions carefully and in their entirety before using the product for the first time. The product must only be used as described in these instructions. Do not use the product for any purpose other than that specified. Please follow all safety precautions and instructions contained herein. The product must only be used by persons who have read and understood these instructions. We assume no liability for damage resulting from improper use, mishandling, or failure to follow these instructions. This does not affect mandatory liability under the Product Liability Act, liability for injury to life, limb, or health, or liability arising from willful misconduct or gross negligence.

Product description

GumQuick Pro is an addition-curing vinyl polysiloxane for the laboratory fabrication of flexible gingival masks on the master model. The material is universally applicable throughout the entire crown, bridge, and implant range. With a final hardness of 72 Shore A, it combines optimal dimensional stability with the necessary flexibility for high-precision work. This avoids non-physiological pressure during wax-up modeling and ensures an optimal fit of the prosthetic restoration. The gingiva-colored design accurately replicates the natural gingival contour and guarantees safe, trouble-free modeling of the secondary parts. The material is fully scannable for the digital CAD/CAM workflow.

Processing

Indirect Method


Apply a matrix of hard, kneadable silicone (e.g., Laborsil pro) onto the dies of the unsawn model, which should already include the gingiva mask. After curing, trim the silicone so the tops of the dies are exposed (approx. 1 mm). The model can then be sawn and the gingiva parts ground as usual. Afterward, verify that the matrix is seated correctly.

To guarantee a flawless surface, apply a separator to the contact areas between the silicone and GumQuick Pro. Apply GumQuick Pro directly from the cartridge onto the model and into the matrix. Position the matrix on the model so that excess material can flow out through the openings at the top of the die. Ideally, the gingiva mask should cure under pressure. After approximately 5:30 minutes, remove the matrix. The mask can then be trimmed using scissors, a scalpel, or rotating instruments.

Direct Method

First, apply DuoSep to the contact surface between the impression material and the GumQuick Pro to ensure a flawless surface.

Caution: To prevent damage to the gingiva mask surface when using GumQuick Pro with polyether impressions, DuoSep must be used! Allow polyether impressions to ventilate for at least 3 hours before producing the mask. Afterward, apply GumQuick Pro directly from the cartridge into the impression around the transfer posts. Thanks to its excellent flow characteristics, the material will securely surround the abutments. After approximately 5:30 minutes, the material is cured, and the impression can be poured with plaster. The gingiva mask can then be easily trimmed with scissors, scalpels, or rotating instruments.

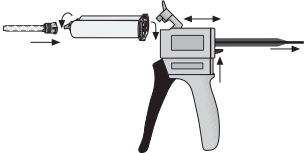Precautions

- Wear protective clothing! Set material is chemically resistant. Stains on your clothing cannot be removed.
- Due to production reasons, components may initially emerge unevenly from new cartridges. For this reason, express a small amount of material before attaching the mixing tip. Also, discard the first increment expressed after attaching the mixing tip.
- Attach a new mixing tip before each use of the cartridge. Using a mixing tip that contains partially or fully cured material may cause the cartridge to rupture, rendering the remaining material unusable.
- Do not use latex gloves, hand creams, cleaning agents, etc. during processing, as the chemical reaction of the two components is inhibited and the material no longer hardens completely. For this reason, we recommend that you wear commercially available polyethylene gloves.
- The setting times can be increased by lower temperatures or reduced by higher temperatures.
- After use, leave the mixing tip on the double cartridge as a cap.
- Use suitable personal protective equipment (respiratory protection, safety glasses, and, if necessary, a suitable exhaust ventilation) during the finishing of the cured material.

Keep out of reach of children!

Disposal

Disposal of contents / container in accordance with local / regional / national / international regulations and safety data sheet.

Hazard identification

Hazard and safety notes can be found on the product label and the corresponding safety data sheet (Download at <https://dentamid.dreve.de/service/downloads>).

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